

Testing Arterial Line

1

Closed loop communication

"Testing arterial line"

2

Centrifugal only → RAP mode on

3A

Pulsatility check

- Remove arterial clamp
- Matches patient pressure

⚠ **Warning** – centrifugal pumps * see CENTRIFUGAL PUMPS

3B

Test Volume

- Give 100 mls FAST
- Observer reservoir volume & line pressure ↑ & ↓ to baseline (MAP)

3C

Reclamp arterial line: Roller + centrifugal

- Centrifugal only, apply to Position 1, engage autoclave, remove Clamp 3


4

Communication

"Line tests well/good"

5

Centrifugal only → RAP mode off

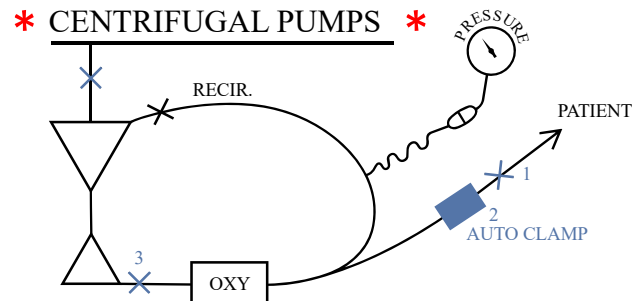

CLAMP SEQUENCE

- ☐ **Add** Clamp 3 (outlet, centrifugal pump)
- ☐ **Remove** Clamp 1
- ☐ **Remove / open** Autoclave #2