

Testing Arterial Line

Verify arterial cannula position and circuit integrity before initiating cardiopulmonary bypass

1

Closed loop communication

"Testing arterial line"

2A

Pulsatility check

Warning
For centrifugal pumps,
see CIRCUIT REFERENCE

- Remove arterial clamp
- Matches patient pressure

2B

Test Volume

- Give 100 mls FAST
- Observe reservoir volume & line pressure $\uparrow$ & $\downarrow$ to baseline (MAP)

2C

Reclamp arterial line: Roller + centrifugal


- Centrifugal only, apply to Position 1, engage autoclamp, remove Clamp 3

4

Communication

"Line tests well/good"

CIRCUIT REFERENCE

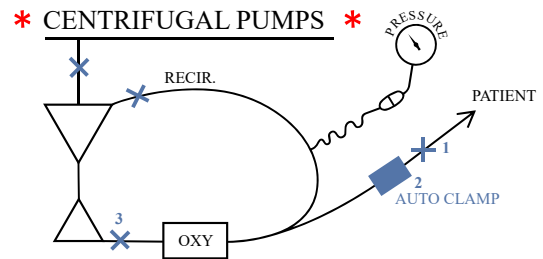

- ☐ Centrifugal only $\rightarrow$ RAP mode on
- ☐ Add Clamp 3 (outlet, centrifugal pump)
- ☐ Remove Clamp 1
- ☐ Remove / open Autoclamp #2

TEST VOLUME - EXPECTED RESPONSE

- ✓ Reservoir volume decreases
- ✓ Arterial line pressure rises
- ✓ Arterial line pressure returns to the patient's MAP
- ✓ No unexpected, sustained elevation in arterial line pressure