

Testing Arterial Line

1 Closed loop communication

"Testing arterial line"

2 Centrifugal only → RAP mode on

3A Pulsatility check

- Remove arterial clamp
- Matches patient pressure

⚠ **Warning** – centrifugal pumps * see CENTRIFUGAL PUMPS

3B Test Volume

- Give 100 mls FAST
- Observe reservoir volume & line pressure ↑ & ↓ to baseline (MAP)


3C Reclamp arterial line: Roller + centrifugal

- Centrifugal only, apply to Position 1, engage autoclave, remove Clamp 3

4 Communication

"Line tests well/good"

5 Centrifugal only → RAP mode off

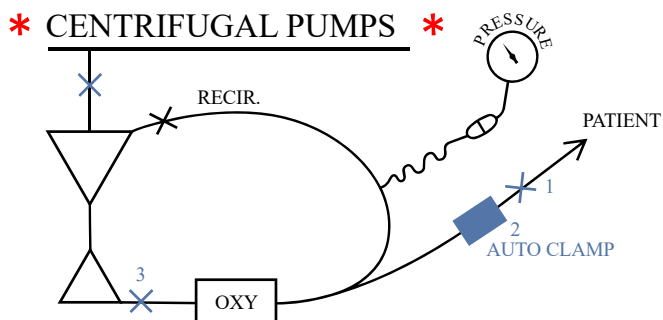

CLAMP SEQUENCE

- ☐ Add Clamp 3 (outlet, centrifugal pump)
- ☐ Remove Clamp 1
- ☐ Remove / open Autoclamp #2