

Testing Arterial Line

Verify arterial cannula position and circuit integrity before initiating cardiopulmonary bypass



1 Closed loop communication


2a Pulsatility check

⚠ WARNING

For centrifugal pumps, see Circuit Reference →

- Remove arterial clamp
- Matches patient pressure

CPB - 2024

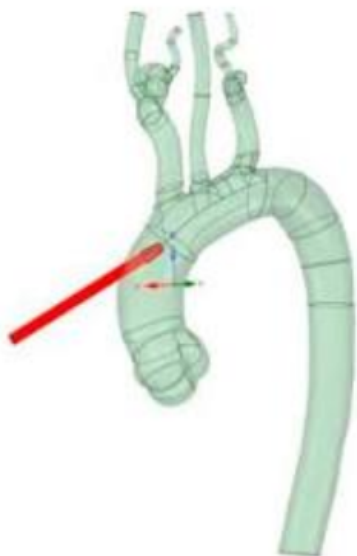

2b Test Volume

- Give 100 mls FAST
- Observe reservoir volume & line pressure ↑ & ↓ to baseline (MAP)

2c Reclamp Arterial Line


Roller + centrifugal

- Centrifugal only → apply to Position 1, engage autoclave, remove Clamp 3

4 Communication

"Testing arterial line"

CIRCUIT REFERENCE

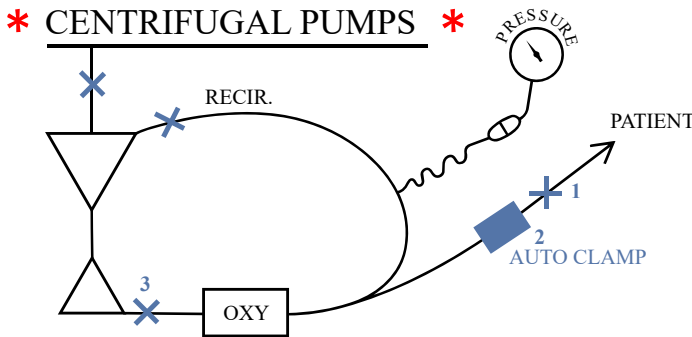

- ☐ **CENTRIFUGAL ONLY** → RAP mode on
- ☐ **Add** Clamp 3 (outlet, centrifugal pump)
- ☐ **Remove** Clamp 1
- ☐ **Remove / open** Autoclave #2

EXPECTED RESPONSE

- ✓ Reservoir volume decreases
- ✓ Arterial line pressure rises
- ✓ Pressure returns to patient's MAP
- ✓ No unexpected, sustained pressure elevation

"Line tests well/good"